

DROPS: CLEANING & LUBRICATION

Periodic lubrication of the drop mechanisms will ensure smooth, crisp drops. Apply **3-IN-ONE** oil every month or 500 adjustments. If the drops have not been lubricated on a regular basis, first apply **WD-40** to clean the plungers from dust & debris. See below.

CERVICAL DROP



Routine lubrication:

1. Elevate the headpiece and cock the drop mechanism.
2. Apply **3-IN-ONE** oil into the hole of the plastic block. Allow a few minutes for the oil to penetrate before using.

If the drop has not been lubricated on a regular basis:

1. First apply **WD-40** into the hole and then drop the headpiece several times to clean the plunger from dust & debris.
2. It is important to follow the cleaning with **3-IN-ONE** oil for lubrication.

TENSION CONTROL ADJUSTMENT

The resistance on the tension control rod can be increased or decreased. Do this by adjusting the set screw on the bottom of the thoracic, pelvic and caudal drop blocks using a 1/8" allen wrench:

Increase resistance: Turn the set screw clockwise. This is helpful if the tension control rod is "backing out" while using the drop.

Decrease resistance: Turn the set screw counter-clockwise.

THORACIC, PELVIC & CAUDAL DROPS



Routine lubrication:

1. Place the table on its cushions.
2. Locate the stainless steel plungers.
3. Apply **3-IN-ONE** oil onto each plunger. Allow a few minutes for the oil to penetrate before setting the table upright.

If the drops have not been lubricated on a regular basis:

1. Lay the table on its cushions.
2. Apply **WD-40** onto the plungers and allow a few minutes for penetration. Set the table upright and drop each section several times to clean from dust & debris.
3. Set the table back on its cushions and apply **3-IN-ONE** oil onto each plunger for lubrication. Allow a few minutes for the oil to penetrate before setting the table upright. Drop each section several times.





Set Screw

HEADPIECE HINGE LOOSENING & LUBRICATION

If the forward motion drop on the headpiece becomes sluggish or only one side is moving, the hinges on each side of the headpiece have likely tightened over time. Loosen as follows:

1. Using a thin 3/8" thin wrench, turn the nut clockwise to the wood, which loosens the hinge.
2. Apply a few drops of 3-IN-ONE oil to the nut & screw on each side of the hinge on a regular basis to keep the hinge lubricated and moving freely.

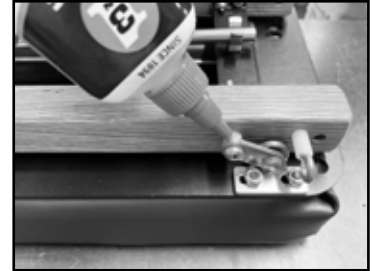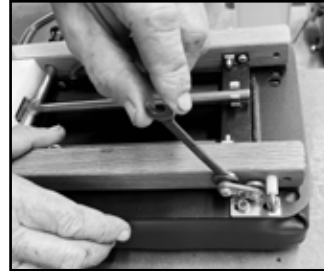

VINYL CLEANING

Water Rinse: It is important to follow any cleaning with a thorough, clear water rinse to minimize premature deterioration of the vinyl from extended exposure to chemicals.

Mild Daily Cleaning:

Use dish soap and warm water (1:10) with a soft cloth. **Follow with a thorough, clear water rinse.**

If more aggressive cleaning is necessary, use a soft bristle brush with the same solution. Avoid harsh detergents and powdered abrasives. Areas coming in contact with hair, body oils or perspiration should be washed frequently. Remove stains immediately to prevent the possibility of becoming permanent.

Disinfecting Options:

Bleach: Dilute 4 tsp bleach in 1 quart of water in a spray bottle. Spray on vinyl, leave for 1-2 minutes and rinse well with water.

Hydrogen Peroxide (3%): Available in spray bottle. Use undiluted and leave on for 5 minutes and rinse well with water.

VINYL REPLACEMENT

In the event of wear or damage, pre-sewn replacement vinyl is available from Thuli Tables.

WARRANTY

Our tables are designed and built to high standards. We are proud to offer a limited lifetime warranty against manufacturing defects, with the exception of the headpiece piston and the motor (on elevation tables), which are covered by a one year warranty.

This warranty is valid to the original owner if the table has not been altered in any way. It does not cover issues caused by normal wear & tear or damage due to accidents, improper use or negligence.

If your table needs repair, whether under warranty or not, please contact us. We will determine what is needed and send you either a replacement or the necessary part for repair. Our tables have been designed so that replacing a part is relatively easy to do yourself. Please note that we do not cover labor costs for replacing a part.